



AUTOMATED CARDIOTHORACIC RATIO CALCULATION AND CARDIOMEGALY DETECTION

[DEEP LEARNING]



Rabia Zafar

[PHDCSF21M503]

College of Information Technology,
University of the Punjab, Lahore

What is Cardiomegaly?

Cardiomegaly is a medical condition in which the heart is enlarged..

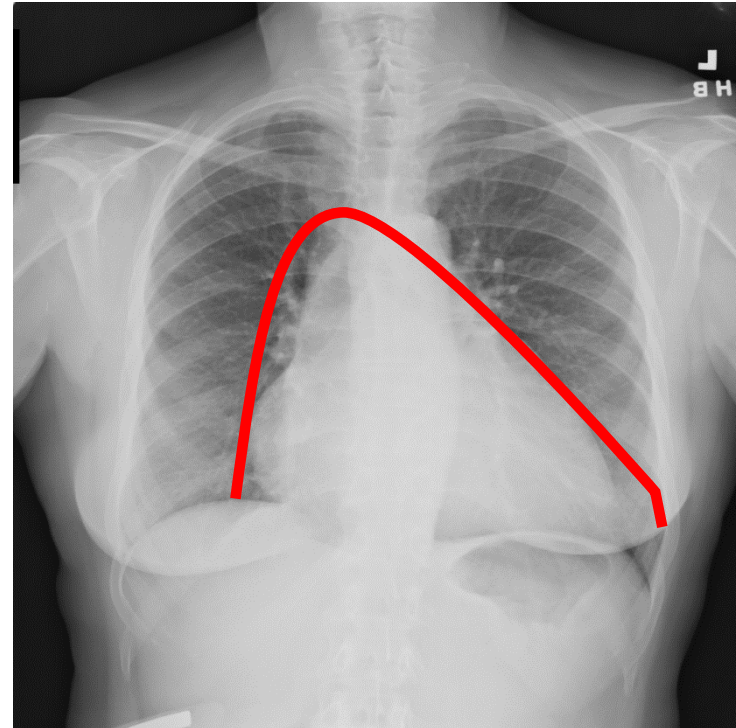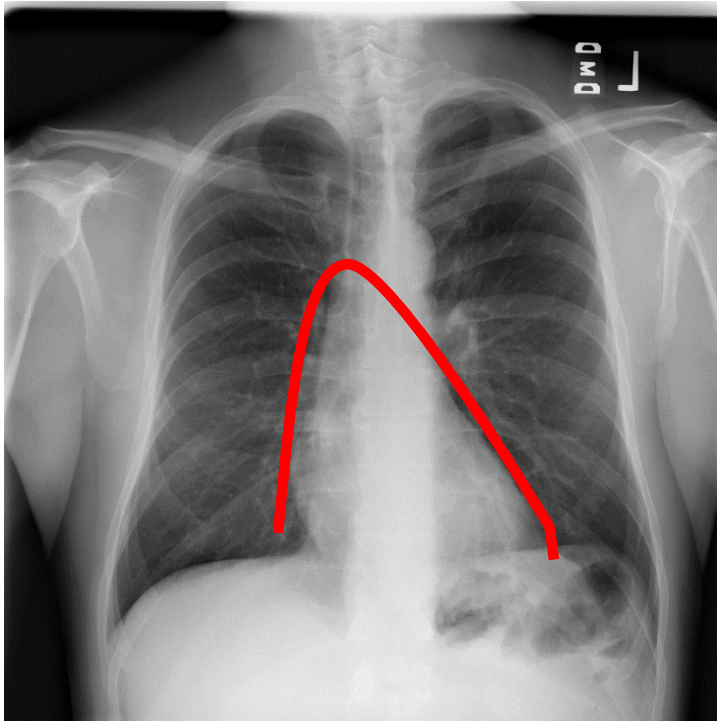

Motivation

- Chest X-ray, or CXR, is widely used in diagnosing abnormal conditions in the chest and nearby structure.
- Radiologists routinely perform cardiothoracic ratio (CTR) measurement
- Picture Archiving and Communication Systems (PACS) include drawing tools to aid the assessment of CTR,
- labor intensive
- time consuming.

CTR

- CTR as one of the most important indicators of cardiomegaly due to the simplicity of the calculation.
- CTR of a chest X-ray image is calculated as cardiac diameter (the diameter of the heart) divided by the thoracic diameter (the diameter of the chest).
- $CTR = (TD)/CD$
- $TD = MRD + MLD$,
- Midline-to-right heart diameter : MRD
- Midline-to-left heart diameter : MLD

Methodology


- Two Approaches are used
- Using Traditional Method
 - Lung Segmentation
 - CTR Calculation
- Using Convolutional Network

Dataset



| Dataset | Total | | Train/Val/Test |
|-------------------------------|--------|------|----------------|
| Kaggle | 342 | | 243/59/40 |
| NIH from Kaggle | 5547 | | 4438/493/616 |
| NIH CHEST X-ray ¹⁴ | ALL | 5550 | 4006/772/772 |
| | Cardio | 1526 | 1006/260/260 |

Lung Segmentation

- Lung Segmentation is done by using pre trained model that uses ResNet34 as Neural Network

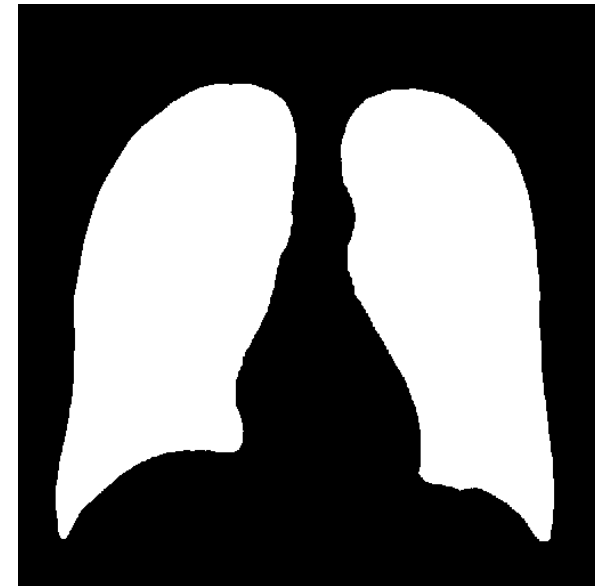

Result of Segmentation

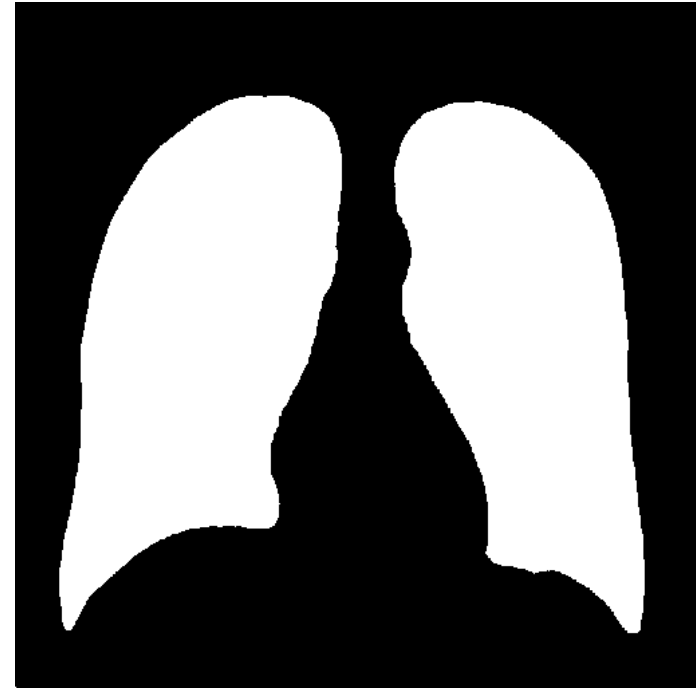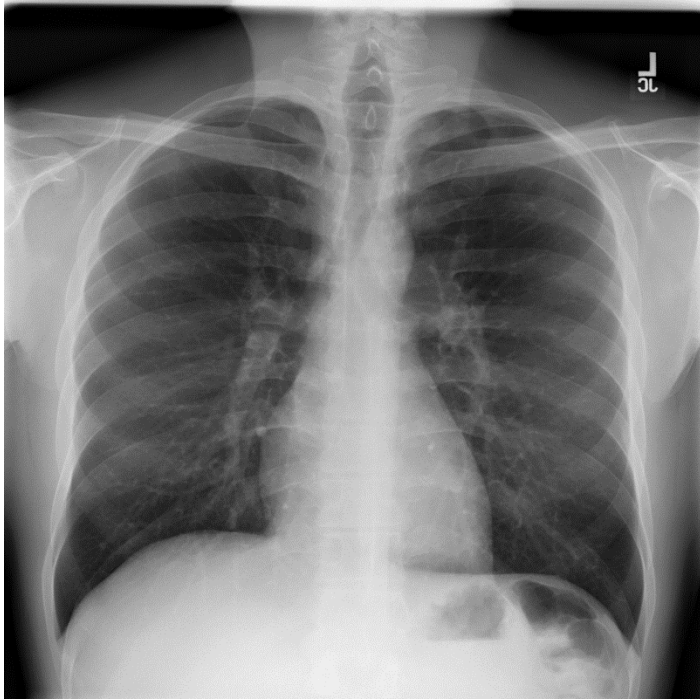

CTR Calculation

- It is done by using traditional image Processing.

Algorithm: CTR calculation

Input-Segmented Image

Output-CTR

Step1: load segmented image

Step 2: Draw a rectangle bounding box around lungs automatically

Step 3: save width of bounding box as CR.


Step 4: Calculate mid in y-direction

Step 5: Calculate $\frac{3}{4}$ in x-direction

Step 6: Calculate MRD and MLD

Step 7: Calculate $CTR = TR / CR$

Methodology

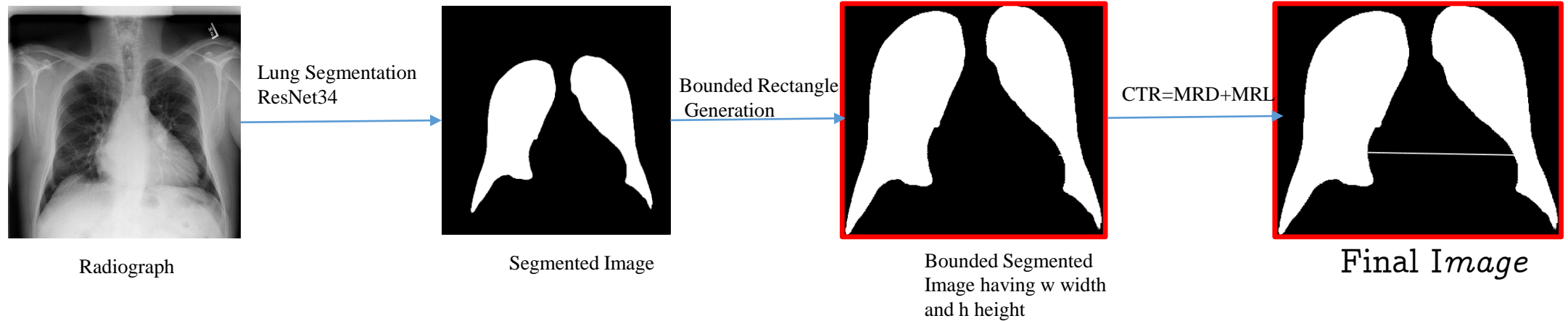

Methodology

- Convolutional Network
 - 3 convolution layer
 - 3 Max-pooling layer (2,2)
 - 5 layer MLP
 - RELU activation function
 - Sigmoid

Results

- Accuracy Attained 70%
- Reason:



Results

| Dataset | Method | Accuracy | | |
|-----------------|----------------------------|----------|------|------|
| | | Train | Val | Test |
| Kaggle | CNN without regularization | 0.10 | 0.91 | 0.80 |
| | CNN with regularization | 0.98 | 0.89 | 0.82 |
| NIH from kaggle | CNN without regularization | 0.99 | 0.73 | 0.74 |
| | CNN with regularization | 0.98 | 0.71 | 0.70 |
| NIH | CNN without regularization | 0.97 | 0.81 | 0.82 |
| | CNN with regularization | 0.86 | 0.81 | 0.78 |

Git-Hub Link

- [PROJECTCV21/Deep-Learning-Project: AUTOMATED CARDIOTHORACIC RATIO CALCULATION AND CARDIOMEGALY DETECTION USING IMAGE PROCESSING AND DEEP LEARNING \(github.com\)](#)

Thank You!